

Operating Manual

3-part Multi Couch on Columns



ISO 13485
Management System Certification
BUREAU VERITAS
Certification Denmark A/S



Operating Manual

3-part Multi Couch on Columns

Content

1. Introduction.....	3
2. Unpacking	3
3. General Security Rules.....	4
4. Area of Use	5
5. Products Life expectancy	5
6. Equipment and Accessories	6
7. Power connection	7
8. Battery charging (accessories)	7
9. Electrical operating procedure.....	8
10. Manual operating procedure.....	9
11. Trendelenburg	9
12. Technical specifications.....	10
13. Handling of parts during Transport	10
14. Cleaning	10
15. Maintenance.....	11
16. Service	11
17. Transport and storage	12
18. Disposal.....	12
19. Electromagnetic – EMC effects	12
20. Guarantee	13
21. Symbol explanation	14
Annex A Periodic inspection.....	15
Annex B Daily inspection before use.....	16
Annex C Recommended annual service inspection	17
Annex D Drawings.....	19
Annex E Spare parts list.....	20
Annex F Troubleshooting	21
Appendix Battery	22

Operating Manual

3-part Multi Couch on Columns

1. Introduction

Congratulations with your new Rehab Care DK ApS equipment. Rehab Care DK ApS products are modern and up to date, based on our long experience in this field. This equipment brings functionality, quality and easy handling to an optimal level. Through continuous development at all levels of our manufacturing process, we are committed to pay attention to details and quality. We are therefore your guarantee of a very reliable and functional product.

The Couch complies with the MDR, classes 1 products (Medical Devices Regulation (EU) 2017/745 - MDR), EN60601-1, EN60601-1-2, IEC21856:2022.

Thank you for your confidence in our product, and for your cooperation.

We wish you every satisfaction with your new equipment.

Best regards Rehab Care DK ApS.

Version 09-09 2024

Manufactured by:



Rehab Care DK ApS
Avnvej 10
7400 Herning, DK

Agent:

2. Unpacking

Unpack the couch carefully and check all parts and accessories.

Check for possible transport damage.








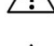





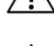
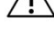
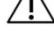

We recommend that the couch is cleaned before use, see cleaning instructions.

Contact dealer if the equipment is damaged or not functioning.

Operating Manual

3-part Multi Couch on Columns

3. General Security Rules

-  Warning! Do not use a high-pressure cleaner on this product.
-  Warning! Do not use chemicals that are not improve. Carry out an annual inspection for damage.
-  Warning! Use only the supplied remote control.
-  Warning! Prevent foreign articles from damaging the wheels.
-  Warning! Use only spare wheels recommended by the supplier.
-  Warning! The couch must not be operated by children under the age of 12.
-  Warning! The couch has a maximum capacity of 250 kg SWL.
-  Warning! Use only approved mattresses of the same type as those supplied.
-  Warning! The couch should not be used if there is any doubt about it being in good working order.
-  Warning! The couch is not approved for private use.
-  Warning! When the patient is not being supervised, the couch should be in the lowest position to prevent falling accidents.
-  Warning! Wrong use of cables, power supply, and remote control can lead to accidents.
-  Warning! Prevent damage to power cable.
-  Warning! Couches must not be stacked on top of one another.
-  Warning! Use of non-approved accessories may result in loss of CE approval.
-  Warning! Attention must be paid when changing the position of the bed guard.
-  Warning! Use of the Trendelenburg function must only be in accordance with the hospital's internal security rules.

Read all instructions carefully before use.

Operating Manual

3-part Multi Couch on Columns

4. Area of Use

Flexible combi gyn./multi-Couch with a variety of adjustments options for treatments.
Used in the hospital and healthcare sector and by private practitioners.

The couch is for indoor use only.

The couches safe working load (SWL) is as follows:

SWL:	250 kg.
Patient:	180 kg.
Mattress:	10 kg.
Accessories:	60 kg.

Persons with physical or psychological handicaps must not use the product without supervision or without adequate training in the use of the couch.

Children must not be left unsupervised around the couch; it is not a toy.



Warning! Couch is not tested and evaluated in relation to use in private homes, therefore, do not use product in private homes



Warning! Couch should not be used by children under the age of 12

5. Products Life expectancy

The product has a life expectancy of 10 years.

Life expectancy is based on the correct use, cleaning, service and maintenance.

If the couch is to be used longer than the above-mentioned life expectancy, it must be evaluated by qualified service personnel before further use.

It is the users' responsibility to evaluate the couch maintains, if used over the listed lift time.

Electrical parts must only be serviced by qualified technicians.

Operating Manual

3-part Multi Couch on Columns

6. Equipment and Accessories

The couch consists of:

- Wheel-based undercarriage with
- Central brakes function on all wheels
- 3-part lying surface
- Electrical remote control
- Electrical height adjustment
- Manual adjustment of back and seat positions
- Trendelenburg/Reverse Trendelenburg adjustments
- Paper roll holder

Extra accessories available:

- Electric adjustment of back and seat position
- Foot control for 1 or 2 functions
- Neck pillow half round
- Cushions in other colors
- Bed guards
- Padded covering for bed guard
- Leg support in plastic cup or padded cup
- Stirrups



Stirrups



Leg support with plastic cup



Leg support with padded cup

Operating Manual

3-part Multi Couch on Columns

7. Power connection

Plug main cable into power supply.

The couch is equipped with a 230V/DC transformer.

Roll up cable when charger is not in use.



Warning! Before connection, control the cable for damage



Warning! Remove cable from main switch to disconnect power supply



8. Battery charging (accessories)

The battery box should be charged before use to ensure the longest possible life of the battery. Charging time depends on the condition of the battery.

Before the couch is used for the first time, the battery should be charged for **24 hours**.

Thereafter battery charging time will be **4 hours**.



Operating Manual




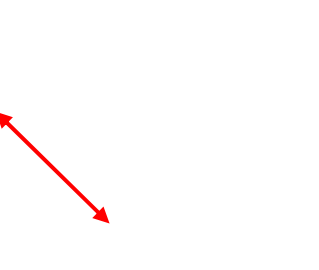
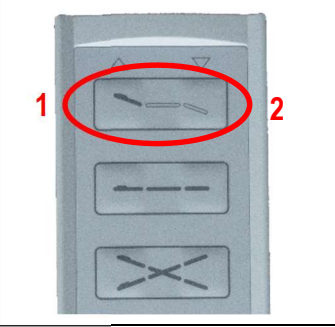
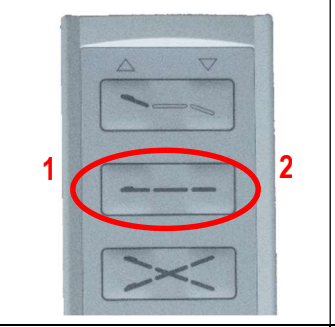
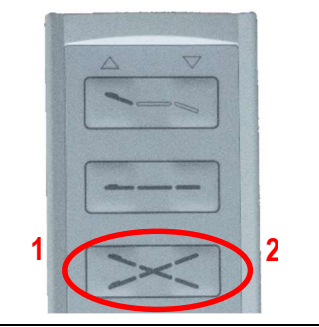
3-part Multi Couch on Columns






9. Electrical operating procedure

Adjustment of the couch is by use of the remote control.

See diagram on the remote control.

The couch can be operated by both caretaker and patient.




			
			
	<p>1. Head part is lifted 2. Head part is lowered</p>	<p>1. Lying surface is elevated 2. Lying surface is lowered</p>	<p>1. Trendelenburg 2. Reverse Trendelenburg</p>

-  When lifting with a patient, the maximum weight must not be exceeded
-  When lifting, SWL must not be exceeded
-  Check there is a free space over the couch before lifting
-  The couch must be able to function unhindered up and down
-  Check the remote-control cable is in good condition and free from any damage

Operating Manual

3-part Multi Couch on Columns

10. Manual operating procedure

<p>Wheel lock (pedal)</p> 	<p>No-lock function</p> <p>Brake pedal in horizontal position – no lock function</p> <p>NOTE: To move the couch, brake pedal should be in horizontal position</p>
<p>Wheel lock (pedal)</p> 	<p>Lock function</p> <p>Brake pedal (red end down) couch is locked against rolling and turning</p>
<p>Wheel steering (pedal)</p> 	<p>Steering activated</p> <p>Brake pedal (green end down) activates steering function (under transport)</p>

11. Trendelenburg

The couch can be adjusted in Trendelenburg (12-degree angle) by the remote control by pressing the button marked with an 'X'.



Warning!
Use of the Trendelenburg function must only be used in accordance with the hospital's internal security regulations

Operating Manual

3-part Multi Couch on Columns

12. Technical specifications

12.1 Choice of mattress:

Madras is attached on the couch.

12.2 Angles and heights:

Angles and heights are given under point **D.1 Couch diagram**.

Lowest	- Measured from floor to mattress top	50 cm
Maximum	- Measured from floor to mattress top	90 cm
Back rest	- Angle	85 degrees
Trendelenburg		12 degrees

13. Handling of parts during Transport

There are no loose parts.

14. Cleaning

The couch should be cleaned with a brush or sponge, using a cleaning or disinfection liquid. Rinse well with clean water and then dry thoroughly.

100% chlorine solution must not be used.

Wheels should only be cleaned with a non-corrosive cleaner.

The mattress should be cleaned with a damp cloth and soap/cleaning liquid.



Warning

Never use a high-pressure cleaner on the product



Warning

Steam cleaning must not be used



Warning

Automatic washing machines must not be used

Operating Manual

3-part Multi Couch on Columns

15. Maintenance

Wheels should be kept free of hair and dirt.

Annual service is recommended.

Moving parts should be checked, tightened and greased.

Worn bearings should be replaced, at turning points.

Batteries should be replaced as necessary.

16. Service

It is recommended that qualified personnel are servicing the couch once a year. It is the responsibility of the owner to ensure the service is carried out. A service protocol for the couch should be established. The protocol should be signed by the customer after each service. If the technician discovers a serious defect with the couch, he is under obligation to remove the couch if the defect could lead to a security risk.

Service, repairs and other work on the couch between service checks should be noted in the service protocol.

16.1 Service - Actuator

The actuator/motor should be replaced after 11.000 cycles.

Replacement of actuator must only be carried out by authorized service personnel.

16.2 Service - Wheels

Wheels should be serviced when necessary: Tightening of bolts.

Dust-tight bearings require no lubrication.

Use only corrosion-free cleaning liquid.

Wheels and parts must be replaced as soon as a malfunction is detected.

16.3 Service - Battery

If the couch is not functioning normally and the batteries are fully charged, it could be that the lifetime of the battery has been exceeded and the battery should therefore be replaced.

Replacement of the battery must be carried out by qualified personnel.

Use only original spare parts.

Operating Manual

3-part Multi Couch on Columns


17. Transport and storage

If the couch is not to be used for a long period, it is recommended that it be covered to protect it from dirt and dust. The battery should also be removed, as the lifetime of the battery is significantly reduced if fully discharged.

The couch should be stored in a frost-free environment with humidity of not exceeding 80%.

Transport and storage requirements:

- Humidity between 10 and 80%
- Temperature between -30 and 40 deg. C
- Pressure between 700 – 1100 hPa

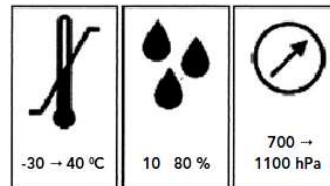

18. Disposal

The couch should be separated into electronic, metal and the mattress waste and should not be disposed to normal waste.

19. Electromagnetic – EMC effects

Mobile RF communication equipment (for example cell phones) can affect electrical medical equipment.

Use of other than the authorized accessories and cables could result in increased emission or reduce product immunity.

Power supply via grid-power must be as normally used in hospital environment.

Operating Manual

3-part Multi Couch on Columns

20. Guarantee

In accordance with Danish law, the couch has a 2-year guarantee.

The guarantee becomes null and void in the following circumstances:

1. Overloading the couch.
2. Violent use of couch.
3. Operating error, wrong use of couch.
4. Lack of or wrong maintenance.
5. If the LINAK electrical parts are opened or have been wrongly used.

At Rehab-Care dk ApS we endeavor to continually improve our products.

We reserve the right to improve the product without prior notice.

All measurements on drawings or other material are only guidelines with the usual provisions for errors.

Operating Manual

3-part Multi Couch on Columns

21. Symbol explanation

The following symbols are used on labels or in this Manual:



Warning, see attached document.



User information



Class 2 equipment



Type B equipment (EN 60601-1)




For indoor use



Medical Equipment.



CE approval MDR (Medical Device Regulation) 2017/745

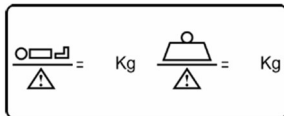

Safe Working Load (SWL) = Patient Max. Weight. + equipment



production address.



Production date.



Product item no.



Product serial no.



unique equipment identifier.

Product label. (example).

The example product label includes the following information and symbols:

- Logo: Rehab Care (RC)
- Production address: Rehab Care, Avnvej 10, 7400 Herning, Danmark
- Production date: 29-07-2024
- Product item no.: REF AKUT-14000
- MAX LOAD: 250 kg
- POWER: 230V/24V5Amp
- DUTY CYCLE: 10 %
- Symbols: Class 2 equipment, CE, MD, Warning, Type B equipment, User information
- QR code
- UDI: UDI
- SN: SN 10013
- MANUAL
- UDI number: (01)05745000344130(11)240729(10)1005(21)10013

Operating Manual

3-part Multi Couch on Columns

Annex A Periodic inspection

A.1 Periodic inspections should be carried out according to the time frame recommended by the manufacturer, or at least once a year. The periodic inspection comprises a visual inspection, especially of the wheelbase structure, hoist mechanism, brakes, turning points, security devices and patient aids, as well as a function test and normal maintenance, e.g. brake adjustment, tightening of screws, etc.

Each inspection must include one (1) lift with maximum load.

A.2 Periodic inspection must be carried out by a suitably qualified person with knowledge of the construction, use and maintenance of the equipment.

A.3 All-important security observations must be entered in a logbook to be kept by the person(s) responsible for servicing/maintenance of the couch.

The date measures are taken in reaction to the observations should also be noted in the logbook.

A.4 Date and result of the inspection should be entered in the logbook and signed by the inspector.

A.5 In case of fault, wear or other damage that threatens the security of the couch, the owner should be informed immediately. In the case of evident security risk, the couch should be withdrawn from use immediately. The couch must not be used until the fault is rectified.

A.6 Any damage threatening security of the couch that occurs between inspections, and which has already been corrected should be entered in the logbook.

A.7 Faults and damage should be reported to the manufacturer/dealer. This feedback must be entered in the logbook. Annex A

Operating Manual

3-part Multi Couch on Columns

Annex B Daily inspection before use

Nr.	Check	Remarks
1	<ul style="list-style-type: none">All parts in place, no loose parts	
2	<ul style="list-style-type: none">No electrical defects, actuator, cables	
3	<ul style="list-style-type: none">No abnormal noise from actuator or control box	
4	<ul style="list-style-type: none">No rust or surface damage	
5	<ul style="list-style-type: none">No visual or other damage	
6	<ul style="list-style-type: none">Wheel and brake function	
7	<ul style="list-style-type: none">Check bolts and mounted parts are properly fixed on couch as well as accessories, in order to prevent undue slackness and spaces	

Daily cleaning

Remove hair and dust around the wheels.

Clean couch with a damp cloth.

Clarification:

The couch should be cleaned with cleaning or disinfection liquid, using a brush or sponge, and thereafter washed. Dry thoroughly after wash.

Concentrated disinfection liquid must not be used.

100% chlorine solution must not be used.

Automatic wash machines must not be used.

For cleaning of wheels, corrosion-free liquid must be used.

The mattress should be cleaned on both sides with a damp cloth and soap.



Warning!
High pressure cleaner must not be used



Warning!
Steam cleaning must not be used



Warning!
Machine wash must not be used

Operating Manual

3-part Multi Couch on Columns

Annex C Recommended annual service inspection

Date:		Product ID/serial no.:
Carried out by:		
Service at:		

No.	Check	Remarks
	Visual check of product	
1	- wear	
2	- distortion	
3	- general damage	
4	- rust	
5	- stability no loose parts	
6	- wheels running easily	
7	- wheel lock functioning	
8	- check bed guard for wear	
9	- check mattress for holes, dirt or damage	
10	- check remote control	
	Lift motor/actuator	
11	- raise couch to maximum to ensure it stops	
12	- actuator/lift motor to be replaced after 11.000 cycles	
13	Check all functions with normal use (full load, complete cycle)	

Operating Manual

3-part Multi Couch on Columns

	Electrical parts/connections	
14	Check all cables for wear	
15	Check all cables to ensure they will not be damaged during use	
16	Check all cables are free of damage	
17	Check all electrical parts are equipped with readable labels.	
18	Check battery and control boxes for wear and tear	
	Other	
19	Clean and lubricate as necessary	
20	Label with date of next inspection	
21		
22		

Date: _____ Initials/signature: _____

The yearly inspection can be electronically registered in the service logbook at Rehab-Care dk ApS.


Contact Rehab-Care dk ApS for further information.

Operating Manual

3-part Multi Couch on Columns

Annex D Drawings

D1. Drawing – Couch

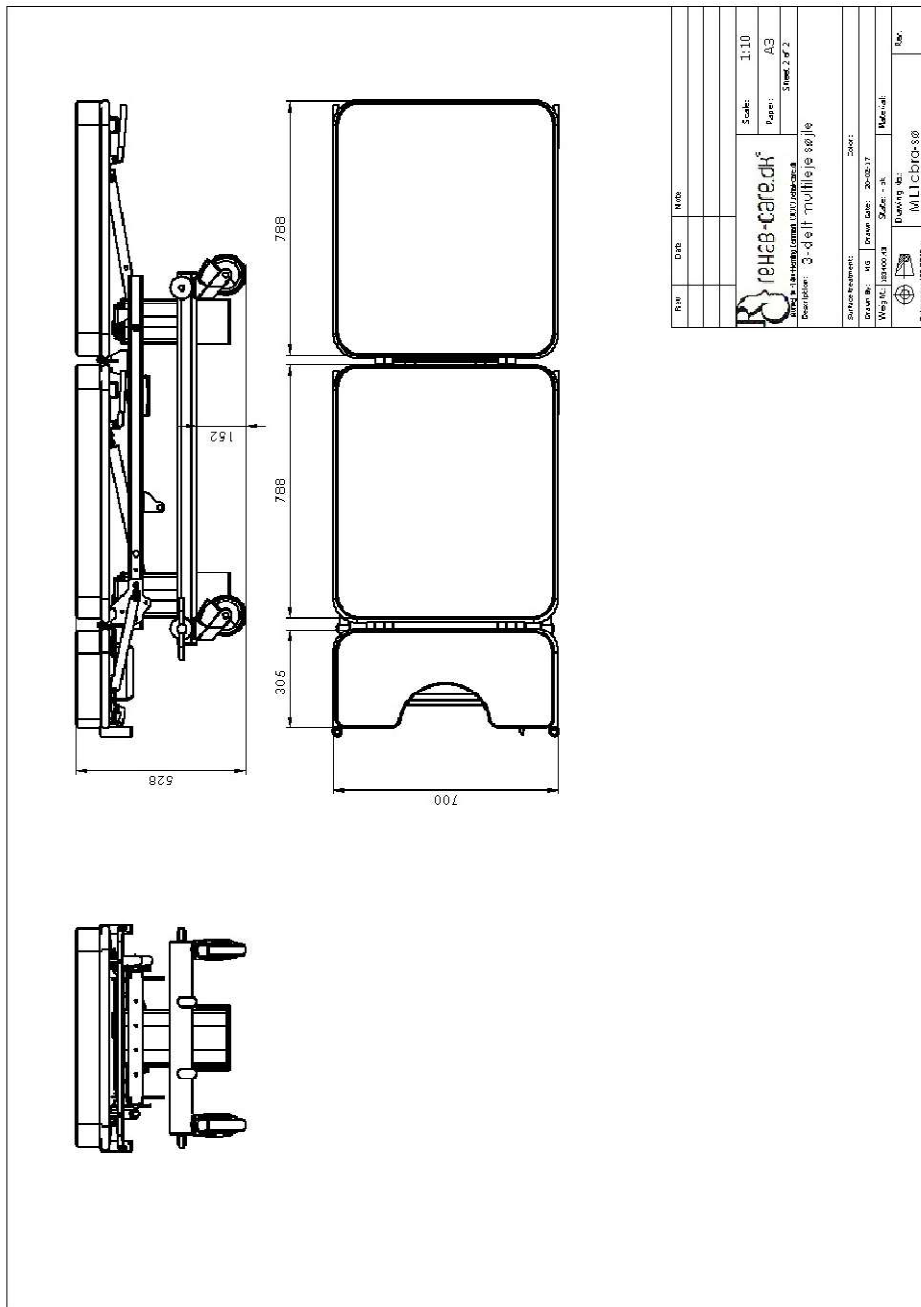


Operating Manual


3-part Multi Couch on Columns

Annex E Spare parts list

ITEM NO.	PART NUMBER	Description	QTY.
1	CG4-100	Gyn sæde	1
2	U5	Rygstykke	2
3	K012RE3-150-195-		2
4	001-500N	Låseshus rullestole	3
5	Forsæd 30x20x6		1
6	L56	Foliesænger til gasstøder	2
7	L56	Udtræksbånd til gasstøder	2
8	GA3-17		1
9	GA3-02		1
10	GA3-02		1
11	WG-200	Lillexsøjle 2 pcs	2
12	L53-H	Ligeflade beslag	2
13	8-1-400-V	Søjle motor MØ1-61	2
14	8-1-400	Søjle motor MØ1-61	2
15	WG-1000	Topramme	1
16	SØ3	Top plade til søjle	2
17	SØ3-1	Aksel for klods	2
18	DIN 912 M10 x 40 ---		1
19	DIN 1537 - M10 ---		1
20	L-4000	Bundramme søjle	1
21	5444-48FJ-125R36-		1
22	L40	Hjul Ø125	4
23	L6-4	Bøsning plusi	4
24	SKA28-1	A-stanser	2
25	L-3010-01	Medbringer bremse	2
26	L-3000-50	Bremse stang	2
27	L42-10	Hjul holder	2
28	L42-20	Bremsestang L=500	2
29	Pedal-højre	Bremsepedal til pedal	1
30	Pedal-venstre	Bremsepedal	1
31	60x30 -lad oval	Plastprop grå 60x30	2
32	Frap Ø16-1,5		2
33	CG-61	Styrebooks	1

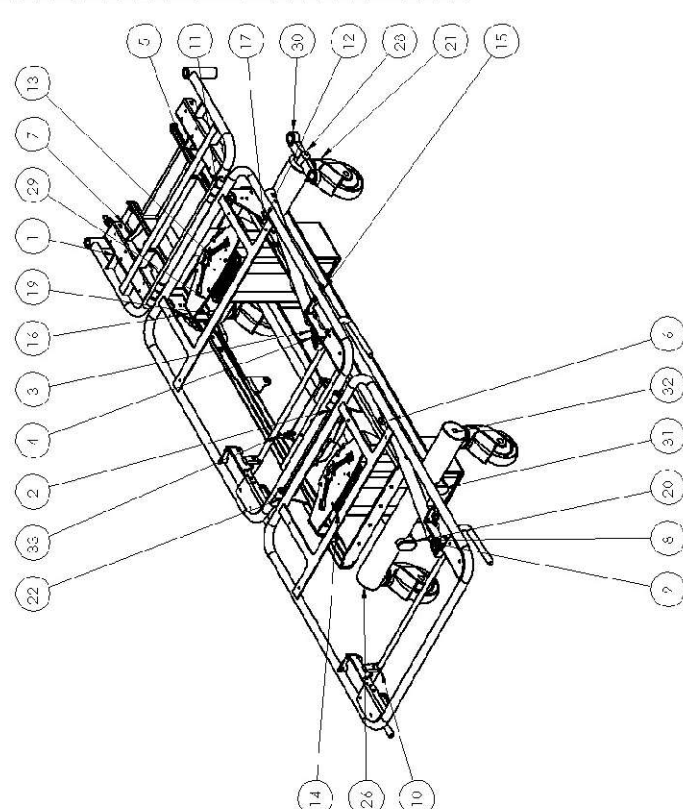

Rev	Date	Desc

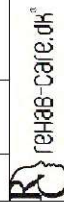

rehab-care.dk
rehabiliteringscenter nordvest a/s

Scale: 1:10
Paper: A3
Sheet 1 of 2

Description: 3-delt multi leje søjle	
Surface treatment: Galv	Color: Galv
Drawn By: MS	Draw Date: 28.02.17
Weight: 8422,87	State: -sk
Material: -sk	Material: -sk
Drawing No: -sk	Material: -sk
Reference: 1902 2718 4L	Material: ML1 CBR SØ

Operating Manual

3-part Multi Couch on Columns

Annex F Troubleshooting

No.	Fault	Possible cause	Solution
1	Malfunction	Is power plug on?	Switch on
2	Malfunction	Are all plugs pressed firmly in?	Press plugs in firmly
3	Malfunction	Is there visible damage to cable?	Replace cable
4	Malfunction	Battery	Replace with fully charged battery

Contact dealer if the couch is not functioning correctly.

A defect couch must not be used.

Operating Manual

3-part Multi Couch on Columns

Appendix Battery

The couch can be equipped with different battery systems.

The following models are used: BA21 and BAT14.

Both systems are approved in accordance with Standard 60601-1 for electrical medical equipment.

Charging

To ensure the longest possible lifetime of the battery box, charge before use.

First-time charging

before the couch is used for the first time, charge the battery for **24 hours** to ensure it is fully charged to ensure the longest possible lifetime for the battery.

Subsequent charging

Subsequently, the battery has a charge time of **4 hours**.

NOTICE:

With storage of longer than 6 months from the production date, the battery must be charged every 6 months.

NOTICE:



The battery can be damaged if it is fully discharged between charges.



Warning! Couch must not be in use during charging



Warning! Battery must not be charged in wet rooms, bathrooms or similar.

BA21	BAT14
	
IPX6	IPX6
Color: grey	Color: White and red
Class 1 electrical safety	Class 1 electrical safety
Use 2 minutes continuously, thereafter 18 minutes disconnected	Use 2 minutes continuously, thereafter 18 minutes disconnected
Working temperature: +5 - + 40 degrees C	Working temperature: +5 - + 40 degrees C